






Arterial Introducer with Haemostatic Valve



CODE	SIZE	COLOUR	DILATOR LENGTH (MM)	DILATOR EXTERNAL DIAMETER (MM)	SHEATH LENGTH (MM)	SHEATH INTERNAL DIAMETER (MM)
453005	5F		155	1,67	100	1,85
453006	6F		155	2,05	100	2,18
453007	7F		155	2,37	100	2,51
453008	8F		155	2,67	100	2,84
453009	9F		155	3,05	100	3,17
453010	10F		155	3,35	100	3,50
453011	11F		155	3,66	100	3,83
453012	12F		155	4,00	100	4,16
453013	13F		155	4,33	100	4,49
453014	14F		155	4,66	100	4,85

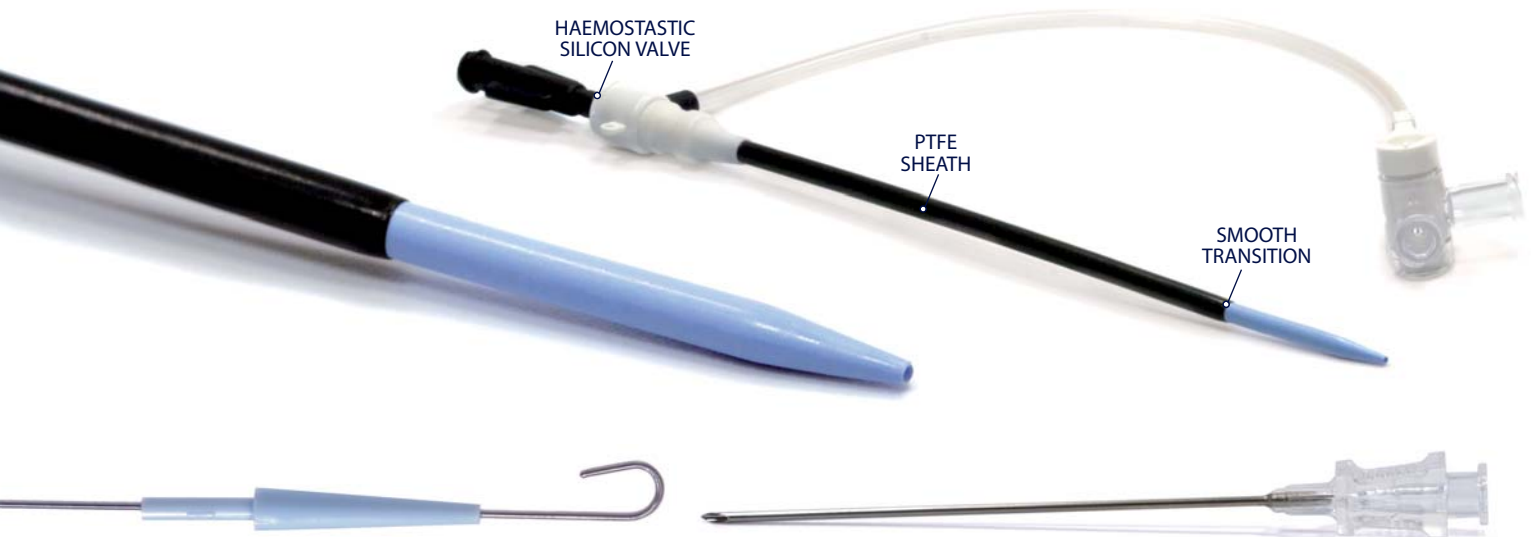

ISOMed Arterial Introducers allow smooth catheter insertion and manipulation. They are equipped with a high performance silicone haemostatic valve that prevents blood reflux and air aspiration.

The dilatator has a long rounded tip for atraumatic insertion. The transition from the dilatator to the sheath is precise and smooth, allowing for the sheath to be easily introduced.

The dilatator locks into the valve assuring insertion without leakage.

Supplied in a sterile (ethylene oxide) blister pack with a J guidewire 0.035" x 45 cm and a puncture needle 18G x 75 mm.

Individual puncture needles, dilatators and sheaths in different diameter and length are available upon request.

Packaging : 10 units/box
Made in France