










Venous Introducer



CODE	SIZE	COLOUR	DILATOR LENGTH (MM)	DILATOR EXTERNAL DIAMETER (MM)	SHEATH LENGTH (MM)	SHEATH INTERNAL DIAMETER (MM)
401004	4 F		140	1.37	85	1.51
401005	5 F		180	1.67	120	1.86
401006	6 F		180	2.05	120	2.17
401007	7 F		180	2.37	120	2.61
401008	8 F		180	2.67	120	2.74
401009	9 F		180	3.05	120	3.21
401010	10 F		180	3.35	120	3.48

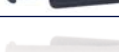

ISOMed Venous Introducers allow smooth catheter insertion and manipulation.

The dilatator has a long rounded tip for atraumatic insertion. The transition from the dilatator to the sheath is precise and smooth, allowing for the sheath to be easily introduced.

Supplied sterile (ethylene oxide), with a straight guidewire 0.035" x 45 cm (except 4F 0.025" x 30 cm) and a puncture needle of 18G x 75 mm.

Packaging : 10 units/box
Made in France